

macafi


Personal

User Guide

Model:
Personal Home



The Healthy Way to Life...

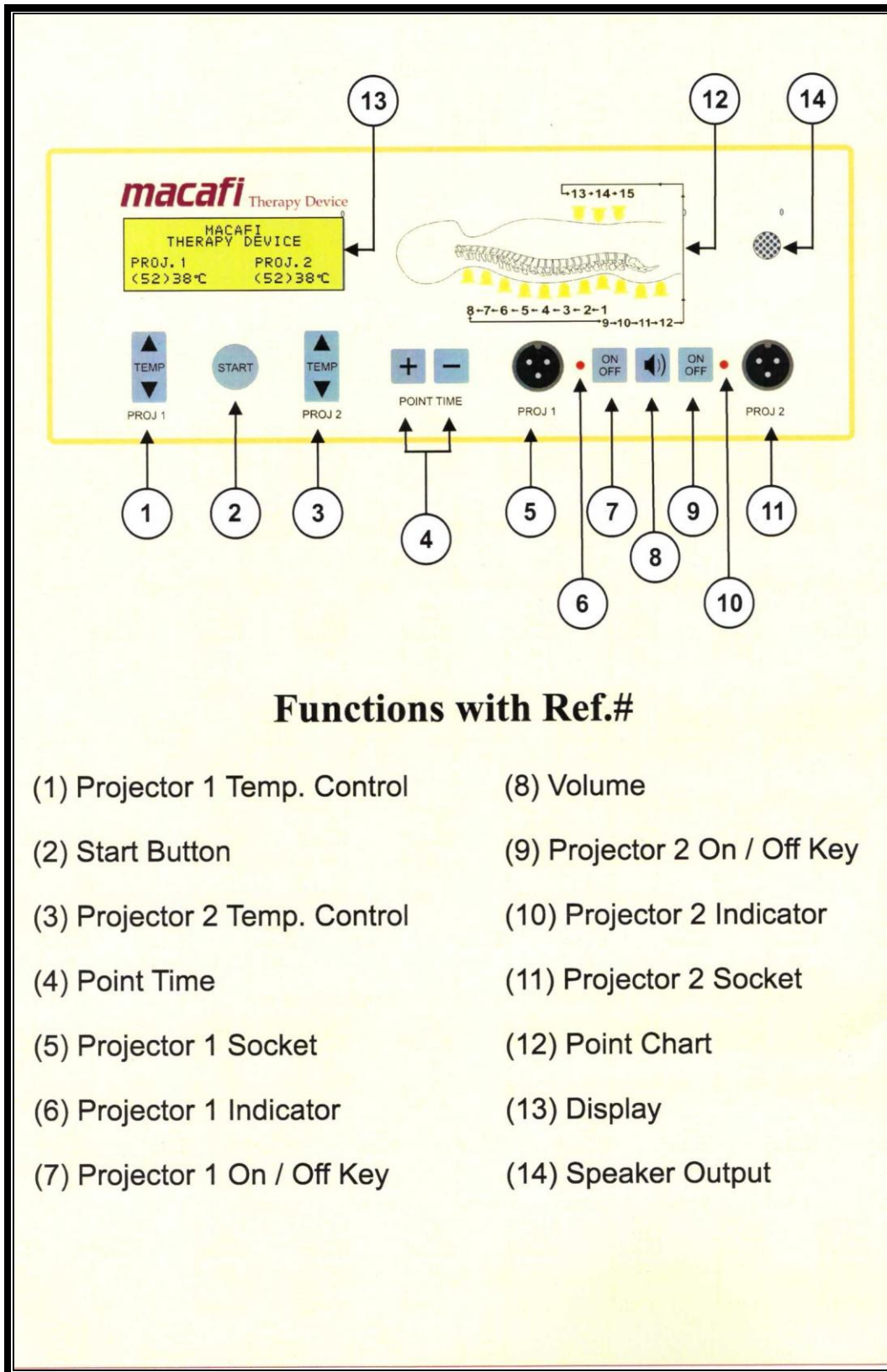

Functions with Ref.#

- | | |
|-------------------------------|------------------------------|
| (1) Projector 1 Temp. Control | (8) Volume |
| (2) Start Button | (9) Projector 2 On / Off Key |
| (3) Projector 2 Temp. Control | (10) Projector 2 Indicator |
| (4) Point Time | (11) Projector 2 Socket |
| (5) Projector 1 Socket | (12) Point Chart |
| (6) Projector 1 Indicator | (13) Display |
| (7) Projector 1 On / Off Key | (14) Speaker Output |

Personal Home



Macafi Controller

Macafi Personal Home is a feature rich product, it comes with one 3 ball and one 9 ball projector, these both projectors can be used simultaneously, this model does not support using two 9 ball projectors simultaneously, two 9 ball projectors are supported only in Personal Ultra model.

Macafi has the feature of voice announcement while therapy is going on, it speaks point numbers at the time of point changing so there is no need to worry about remembering of which point is going on, just relax and get healthy.




9 ball and 3 ball projector

All major operations are done with the front panel keyboard, each and every object of the front panel is shown here on the right page with reference # and its function is described on the following pages with relative reference numbers.

Power On

How To turn ON the Macafi Controller :

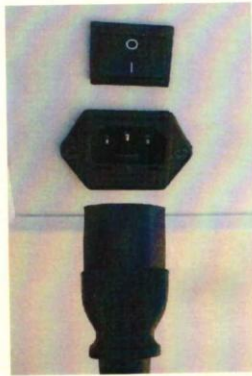

As shown in the picture, there is a main power switch and power cable connector at back panel, plug in the main power cable and turn the main switch on.

Turn off the main power switch when Macafi is not in use.

Turn OFF the main switch at the back to get the Macafi Controller switched off completely.

(1) Projector 1 Temperature Key



Proj. 1

Temperature for Projector 1 can be set from 30 C to 60 C with this Key.

Press the Up arrow key to increment Projector 1 temperature by 1 C.

Press and hold it to continue incrementing

Press the down arrow key to decrement Projector 1 temperature by 1 C, Press and hold it to continue decrementing.

Set temperature of Projector is shown in the bracket and actual temperature next to the set temperature without bracket.

(2) START Points



START Key

Press this key to start the therapy, point counting will start and the current point number and point time will be shown on the display, the point numbers will be announced by voice depending on the volume status (Ref. #8)

Pressing this key while the points are going on will stop the point and reset the point count.

(3) Projector 2 Temperature Key



Proj. 2

Temperature for Projector 2 can be set from 30 C to 60 C with this Key.

Press the Up arrow key to increment Projector 2 temperature by 1 C.


Press and hold it to continue incrementing

Press the down arrow key to decrement Projector 1 temperature by 1 C, Press and hold it to continue decrementing.

Set temperature of Projector is shown in the bracket and actual temperature next to the set temperature without bracket.

(4) Point Time

Point Time

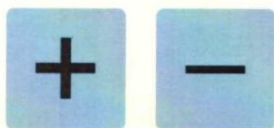

Point Time Keys

Press + key to increment Point Time by 5 seconds.

Press and hold it to continue Incrementing.


Press - key to decrement Point Time by 5 seconds.

Press and hold it to continue decrementing.

Point Time is a value of time duration to be taken for each acupressure point, The time duration can be set from 20 seconds to 180 seconds.

(5) Projector 1 Socket


Projector 1 Socket is connected to the cord of Projector 1.

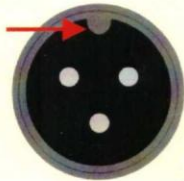

How to connect:

Match the half round tongue on top of the socket (shown by Red Arrow) with the slot in projector connector, insert the connector firmly to this socket and fasten it by turning the ring of projector connector.

(6) Projector 1 Indicator

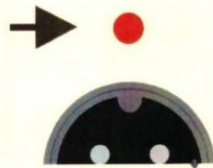

Proj. 1 Indicator

This Indicator lamp shows the status of Projector 1 , it remains ON during the normal operation when Projector 1 is enabled, it goes off when Projector 1 is manually switched off by Projector 1 On / Off key (Ref.#7)

Or when the projector is switched off by the system during the projector error.

(7) Projector 1 On / Off

Projector 1 can be switched off manually by pressing this key.



Projector one lamp will turn off when the projector is manually switched off.

Projector 1 can be switched on again by pressing this key, it will turn on automatically when start key is pressed to start the point.

(8) Volume



Volume Key

Press this Volume key to change the sound level,

Each press of this key will turn the volume in a sequence of Hi, Low and Off, the current status will be announced by the voice.

(9) Projector 2 On / Off


Projector 2 can be switched off manually by pressing this key.



Projector one lamp will turn off when the projector is manually switched off.

Projector 2 can be switched on again by pressing this key, it will turn on automatically when start key is pressed to start the point.

(10) Projector 2 Indicator

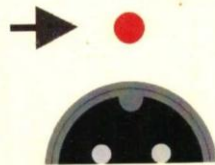


Proj. 2 Indicator

This Indicator lamp shows the status of Projector 2, it remains ON during the normal operation when Projector 2 is enabled, it goes off when Projector 2 is manually switched off by Projector 2 On / Off key (Ref.#9)

Or when the projector is switched off by the system during the projector error.

(11) Projector 2 Socket


Projector 2 Socket is connected to the cord of Projector 2.

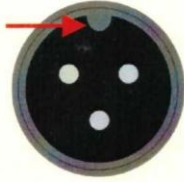


How to connect:

Match the half round tongue on top of the socket (shown by Red Arrow) with the slot in projector connector, insert the connector firmly to this socket and fasten it by turning the ring of projector connector.

(12) Point Chart

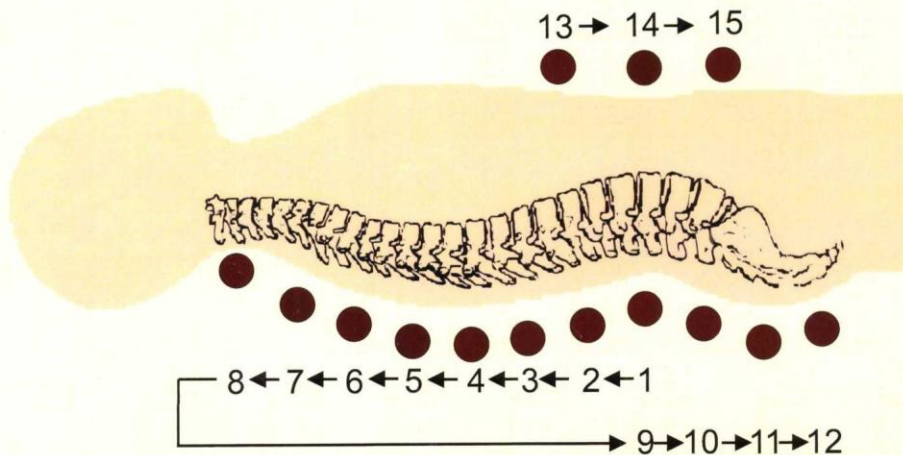


(13) Display

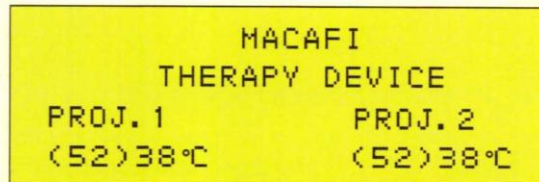

The Liquid Crystal Display is an alphanumeric display that makes the Macafi Controller more user friendly by showing the operational menu and other settings, point time, number of points, projector temperature etc. in a simple way which is very easy to understand.

(14) Speaker Output

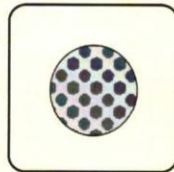

This is a speaker window. When we start therapy, point announcement sound will come out from this window.

Do's and Don'ts

Disconnect the power plug from the wall outlet if the product is not in use for a long time.

To avoid electric shock or fire, take care not to bend the electric cables or wires, or place them underneath heavy sharp objects.

Do not plug too many electric appliances into one outlet as it may cause spark or fire.

Do not place any combustible object (such as a lit candle, cigarette or alcoholic drink), nor any heavy objects or water on top of the product as it may cause electric shock, fire or damage to the product.

To avoid electric shock or fire, do not dismantle, repair or alter the product in any way.

Stop using the product in the following circumstances and contact the nearest center or company to check if:

- Any cord / plug or connector is damaged.
- The product is physically damaged.
- The product is wet.

Do's and Don'ts

Consult your doctor before using Macafi if you are using a pacemaker or any other life supporting apparatus or if you are suffering from fever, skin problems or if you have gone through a medical surgery or having other abnormality.

Do not use Macafi during the pregnancy.

Use cotton clothes and adjust the temperature to suitable level before getting started, you can vary the temperature within the range depending on your body condition and the thickness of your clothes. improper temperature may burn the skin.

This product has a heated surface, persons insensitive to heat must be careful when using the product.

Children should not use the product without adult supervision.

Do not keep the projector on one point for more than 3 minutes, however the same point can be repeated after taking the other points.

Do not take bath immediately after the therapy

It is better not to get the therapy immediately after the meal.

Do not touch the controller or projector or power plug with wet body or wet hands.

Trouble Shooting

Projector Does not work:

Check if it is manually switched off. For more detail, see Ref. #7 for Projector 1 and Ref. #9 for Projector 2.

Macafi Personal Home Model has two Projector sockets, one for 3 Ball projector (Ref. #5) and one for 9 Ball Projector (Ref. #11), If 9 Ball Projector connected to 3 Ball Projector Socket, it will turn off and show the error message on display - "Invalid Projector"

Connect it to relevant 9 Ball socket to get it working, you may have to manually turn on to remove the error message, See Ref. #7 and #9 for more details on how to.

If display shows- "3 Ball Proj.Error" or 9 Ball Proj. Error"
Switch off the Macafi Controller and switch on again, if problem persist again, consult the nearest Macafi Center or company.

No Sound announcement:

Check if volume is set by pressing the volume key,
See Ref. #8 for more details.

If problem persist, Please consult the nearest Macafi Center or Company.

Unable to get the Power On:

Check if the power cable is properly connected at back panel.

Check if the main power switch is on at back panel.

If problem persist, please consult the nearest Macafi Center or Company.